91*181mm

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
250213

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XTherm

**TERMOMETRO
MULTIFUNZIONE 3in1**

**THERMOMETER
3in1 multifunction**

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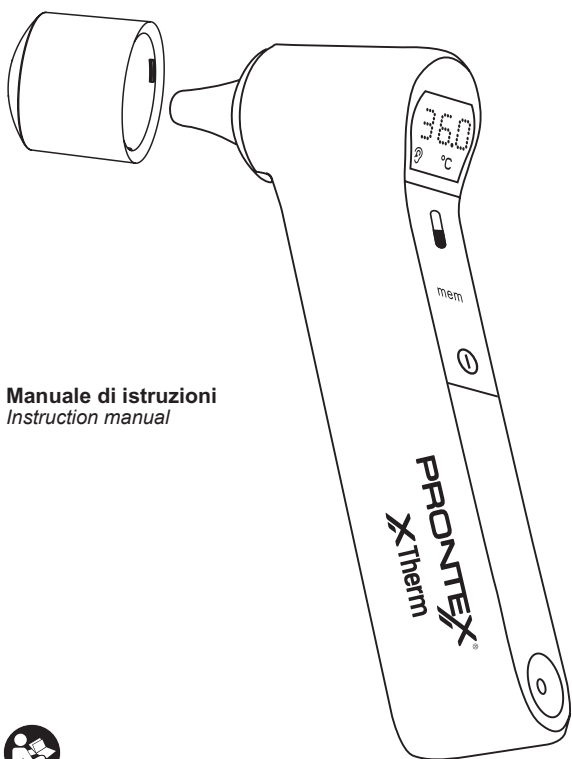CATALOGUE

PRECAUTIONS BEFORE USE	16
DEVICE DESCRIPTION	18
USE CAUTION ITEMS	19
COMMON SENSE ABOUT BODY TEMPERATURE	20
PRODUCT LAYOUT	21
THE SCREEN DISPLAY MARKS DESCRIPTION	21
BATTERY INSTALLATION EXPLANATION	22
BASIC PARAMETER INSTRUCTIONS	22
INTRODUCTION OF MEASUREMENT METHODS	23
PRODUCT CLEANING METHODS INSTRUCTIONS	25
FREQUENTLY ASKED QUESTIONS AND SOLUTIONS	25
TROUBLESHOOTING	26
SPECIFICATIONS OF THE PRODUCT	26
WARRANTY	29

Precautions before use

- The warning signs and illustrations indicated in the manual are intended to enable you to use the product safely and correctly.
- Warning signs, illustrations and their meanings are as follows.

Caution: please refer to attached file.

Legends	
	This mark means warning things (the things that you have to follow), the left graph shows general compulsion.
	Caution
	This mark means prohibited things (not allowed things), it means general prohibition in left picture.
	This mark means prohibited disassembly.
	Type BF Applied part
	It means that you need to read the manual carefully before use.
	It means the package of this unit can comply with the requirements of green environmental protection.
	It means that the material of the device or the device itself is made of renewable material so that we can cycle this product.
	No discard it at will is allowed.
IP	Classification: IP22
	Stacking limit by number
	This way up
	Manufacturer
	Importer
	Indicates the Authorized representative in the European Community.
	Symbol for CE Mark. This symbol certifies that a product has met European Union consumer safety, health, or environmental requirements.

CAUTION

- The patient must follow the doctor's instructions.
- Press the " ⓘ " button and remove the device in case of discomfort.
- Please ask professional doctors to explain the measured value of body temperature.
- The device is intended for measuring human body temperature as well as the surface temperature of objects and the environment
- Please do not use this device for purposes other than body temperature measurement.
- This device contains sensitive electronic components, and should not be used directly in an environment with strong electromagnetic interference (such as near mobile phones, microwave ovens, etc.) to prevent temporary impact on its accuracy.

- Please do not disassemble or repair this device by yourself including changing device. **Any unauthorized intervention will result in the invalidation of the warranty.**
- It is forbidden to collide, drop or mix this device with sharp objects.

**CAUTION**

- This device is only used for human body temperature measurement without use of disease diagnosis; it cannot be used for emergency and continuous measurement in surgery.
- The kids should get away from this device. If the kids are in this range from zero to twelve years old, they need to use this unit to measure under their parents help.
- The patients cannot diagnose the disease and get treatment by themselves on the base of measurement result, they must follow the instruction of doctors.
- Children under the age of 12 and those who are unable to express their thoughts are prohibited from using it on their own.
- This device cannot be used with HF surgical equipment.
- Do not use this device for these people who suffer from otitis external, tympanitis and other ear diseases.
- It is forbidden to immerse the device in any liquid, and it is forbidden to use it in an excessively high or low temperature environment for a long time.

If you use or store this device beyond the range of specified temperature and humidity, It may not ensure accurate measurements.

Use environment:

- Temperature: from +10°C~+40°C
- Humidity: from 15%RH~93%RH

Storage environment:

- Temperature: from -25°C~+55°C
- Humidity: from 0%RH~93%RH

DEVICE INTRODUCTION

Intended use: Prontex X-Therm intended to measure human body temperature by measuring ear canal or forehead. The device can be used on people of all ages at home or in hospital environment, children under 12 years of age and people who are incapable must be assisted.

Scope of application: The device displays body temperature by detecting thermal radiation.

Features:

1. Non-touching type infrared measurement of ear/ forehead temperature.
2. Multiple colors and backlight display: White, Green, Orange and Red.
3. 9 sets of memory values.
4. The switch of degree Fahrenheit °F and degree Celsius °C.
(original setting is degree Celsius °C)
5. Instant measurement in 1 second.
6. Practical and convenient design with no need for probe covers.
7. It has the function of sound on/off.
8. The device idle time of 30 seconds, turn off power automatically.

Warm tip: The measurement results of this unit at any time only can be a reference, it can not replace the medical diagnosis of professional doctor. If you have any question about the personal temperature measurement result, please use and diagnose under the right instructions of doctor.

CONDITIONS FOR ACCURATE MEASUREMENT

CAUTION

1. Consult your doctor to interpret the measurement results. **NEVER change a prescribed medication or its dosage without first consulting your doctor.**
2. Please do not touch with your hands or blow infrared sensor with your mouth. When the infrared sensor is damaged or dirty, it may cause abnormal measurement results.
3. To avoid abnormal readings, allow the device to reach room temperature (measurement location) for approximately 30 minutes before taking the measurement.
 - **Otherwise, may result in incorrect measurement results.**
4. Prevent the child from measuring their temperature on their own.
 - If a child attempts to measure their temperature alone, they may risk causing ear damage.
 - If he swallows the battery or the transparent cover accidentally, please contact your doctor immediately.
5. When measuring body temperature, please do not get close to air conditioning.
 - Avoid affecting the measurement accuracy.
6. Before and after use every time, use a cotton swab dipped 75% alcohol to wipe the probe surface. (If you see stains, fog or water on the infrared sensor glass, please use a cotton swab dipped 75% alcohol for wiping the infrared sensor glass gently.)
 - If you wipe it with disposable paper, it will scratch the infrared sensor resulting in incorrect measurement result.
 - Clean thoroughly to avoid cross-infection of ear disease.
7. The device may suffer mechanical damage.
8. The device touches water or immerses water accidentally, please fully dry before use, especially the water on the surface of the sensor should be clean by using cotton swab.

CAUTION

1. Please do not use this device for these people who suffer from otitis external, tympanitis and other ear diseases
 - It is possible to worsen the pathology.
2. Please do not use this device after swimming or bath or wet ears.
 - It is possible to have low measured temperature value.
3. Dispose of the batteries in accordance with local regulations.
4. If there is allergic reaction after use the device, please stop using it and consult with your physicians.
5. A notice to the user and/or patient that any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

SUGGESTIONS

1. When you tell the doctor measured body temperature value, please tell the doctor that you measured it with ear thermometer.
2. Please don't force to impact, fall, trample and shake this product.
3. Please do not disassemble, repair and modify this device.
4. Please do not allow liquid (such as alcohol, water- drop, hot water and so on) to enter the device body because of this device without water resistance.
5. The product must be kept clean in a dry place.
6. If you find any problems, please contact the sales, you cannot repair the device by yourself.
7. Please do not use it under the environment of electromagnetic interference.
8. Please deal with the waste and residue of this device at the end of the service life according to local laws and regulations.

COMMON SENSE ABOUT BODY TEMPERATURE

The comparison of different measurement methods.

The measured values are different if we use different measurement methods. The WHO provides normal human body temperature reference values, please see below table about the specific temperature difference.

Measurement methods	Normal body temperature
Anal temperature	36.6°C~38°C
Oral temperature	35.5°C~37.5°C
Axillary temperature	34.7°C~37.3°C
Cochlear temperature	35.8°C~38°C

The changes in human body temperature:

Body temperature varies throughout the day.
Below are the details:

At night

Lowest Body temperature is lowest because of sleep and decreasing activity.(below 37°C).

In the morning

Higher From warm bed to the lower temperature room in the morning, the whole body's muscles get contractions and produce heat.

At noon

Highest After lunch, human body reaches the highest temperature and the body will adjust naturally.


Three or four o'clock in the afternoon

Lower Due to physical exertion, blood sugar decreased.


In the evening

Lowest Due to the sun down, room temperature goes down.

PRODUCT LAYOUT

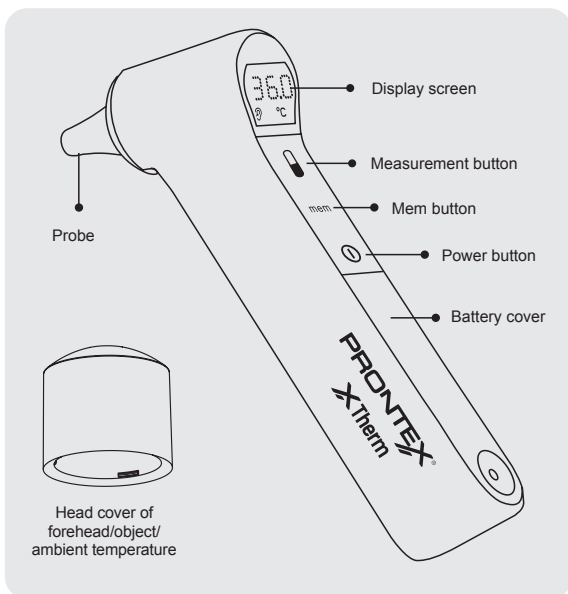
THE SCREEN DISPLAY MARKS DESCRIPTION

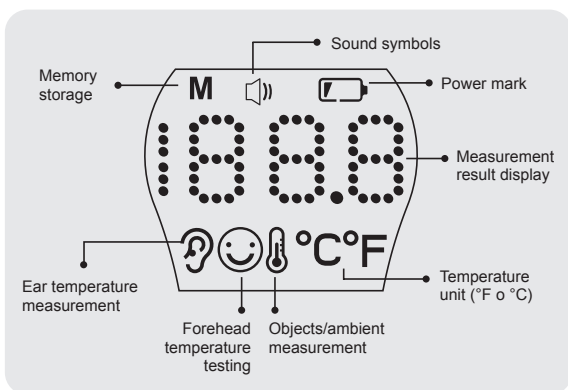BATTERY INSTALLATION EXPLANATION

Battery installation flow:

1. Press battery cover, the battery cover will bounce automatically.
2. Prepare two pieces of battery with the module of 1.5V AAA (number seven battery) batteries (It is recommended to use alkaline batteries), please install it into base of battery according to correct positive and negative poles (as showed inside).

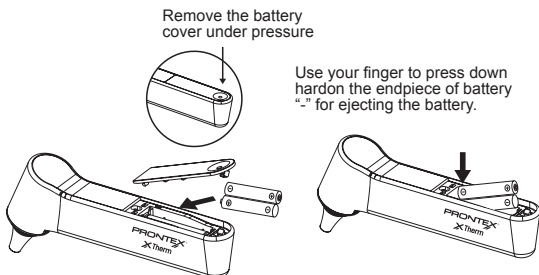Insufficient power caution:

when battery level is low, shows "LO" sign and battery mark is always on, which means you should replace the battery.

Warm prompts

- If you do not use this device for a longtime, please take out battery to assure its longevity. The liquid leakage of battery will harm the device; as well as pollute our environment.
- It is recommended to use alkaline batteries.
- The way to handle scrapped battery should match the requirements of local country government and environmental protection organ.
- **Dispose of used batteries separately from the device.**

BASIC PARAMETER INSTRUCTIONS

1.Sound function: turn on/off

- 1) In the power-on state, press "mem" key to set the sound on or off.
- 2) Press "mem" button, LCD screen will show "🔊", it means the sound works, meanwhile we will hear a short beep sound.
- 3) Press "mem" button again, "🔊" will disappear, it means the sound is off.


2. The switch of °F or °C

In shutdown state, press and hold the "mem" button for 6 seconds, when °F or °C flashing, press mem button again to change unit to °C or °F. Wait for 8 seconds to turn off automatically or press "🔊" to turn off.

MEASUREMENTS STORED IN MEMORY WILL ALSO BE CONVERTED TO THE SELECTED UNIT OF MEASUREMENT.

3. Memory storage function

In off status, press "mem" button, the device can read and save 9 sets of measurement values in order (as below picture shows). It will turn off automatically without operating for 30 seconds or press "I" button to turn off this device with your hands.

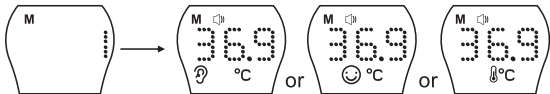4. Back light status instructions

When measured temperature is $<34^{\circ}\text{C}$, it shows LO with red backlight.

When measured temperature value is $34^{\circ}\text{C}\sim 37.1^{\circ}\text{C}$, it shows green backlight.

When measured temperature value is $37.2^{\circ}\text{C}\sim 38.1^{\circ}\text{C}$, it shows orange backlight.

When measured temperature value is $38.2^{\circ}\text{C}\sim 43.0^{\circ}\text{C}$, it shows red backlight.

When measured temperature value is $>43.0^{\circ}\text{C}$, it shows red backlight and shows HI.

Warm prompts: This function is for reference only.

INTRODUCTION OF MEASUREMENT METHODS

Children up to 12 months

Gently pull the child's ear back slightly.

Ear temperature

Children and adults

Pull the ear back and upward.

Ear temperature

Children and adults


Forehead

Forehead temperature

1. Ear temperature measurement

Warm prompts:

Remove the head cover of the device which can measure ear temperature.

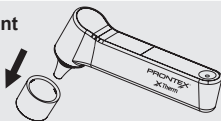1.1 After removing head cover, please press "I" button to turn on this device enter the ear temperature measurement mode, it will display measurement value last time. Then put the thermometer in your ear, the ear temperature can be measured directly by pressing "P" button.

1.2 After 1 second, you can see the measurement result.

When display is flashing, We are waiting for starting measurement

Measurement result


***Fever alert: short beeps for ten times come with urgency.**

Note: If you did not hear the beep sound, which represents the temperature measurement has not yet been completed. Please do not remove the device probe from ear canal at this time. **(If you have closed the sound indication, it will have no sound indication).**

2. Forehead temperature measurement


Warm prompts:

Install head cover of the device which can measure forehead temperature.

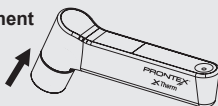2.1 After installing head cover of the device, press “**1**” button to turn on this device to enter the forehead temperature measurement mode, display will show measurement value last time. Then the device is aimed at the forehead and measurement distance should be 0~5mm, the forehead temperature can be measured directly by pressing “**2**” button.

2.2 After 1 second, you can see the measurement result.

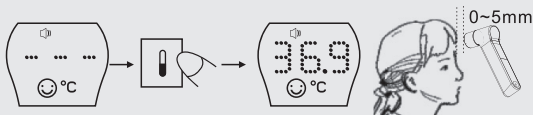When display is flashing, We are waiting for starting measurement.

Measurement result.

*Fever alert: short beeps for ten times come with urgency.

Note: If you did not hear the beep sound, which represents the temperature measurement has not yet been completed. Please do not remove the device probe from forehead at this time. **(If you have closed the sound indication, it will have no sound indication).**

3. Object pattern measurement

Warmprompt:

The installation of the head cover of forehead temperature has no influence to object pattern measurement.

3.1 Long press “**1**” button “**1**” to enter the object measurement mode for 6 seconds, and then aim the thermometer at the object, the temperature of the target object can be measured directly by pressing “**2**” button.


When display is flashing, We are waiting for starting measurement.

Measurement result

Note: If you did not hear the beep sound, which represents the temperature measurement has not yet been completed. Please do not remove device probe from target object at this time. **(If you have closed the sound indication, it will not have sound indication).**

PRODUCT CLEANING METHODS INSTRUCTIONS

In order to ensure the precise measurements, it recommended to clean the device after each use. Please use cotton swab to dip 75% alcohol to wipe the thermometer probe and remove the residue and dirt. We recommend you to clean the probe at least for 5 seconds. Then wipe the other components with a soft and dry cloth for 10 seconds at least. After cleaning, please ensure that there is no visible dirt or residue on your device. Wait 10 minutes before measuring. The device can be cleaned for maximum 40000 times.

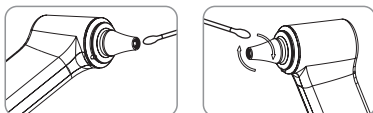Warm prompt:

Why do you need to clean the device after each measurement?
Because infrared temperature adopts highly sensitive technique to detect the temperature of the target object, not only will any earwax and dust paste affect the measurement accuracy, more likely to cause bacteria infection. So we suggest that you should clean your device as the picture shows after you use every time.

FREQUENTLY ASKED QUESTIONS AND SOLUTIONS

Screen display	Reasons	Solutions
HI	When the target object temperature is higher than measurement range which the ear temperature is more than 43.0°C, the screen will show "HI" indication.	When the measurement probe is not properly placed in the ear canal or if the measurement distance is too far during the measurement process, the measurement result may be not accurate.
Lo	When the target object temperature is lower than measurement range which the ear temperature is less than 34°C, the screen will show "Lo" indication.	When measurement probe is dirty, the measurement value may be not accurate; you should use a cotton swab to dip alcohol to wash the measurement probe properly.
Er.H	This device has an upper limit of operating temperature of 40°C. When the ambient temperature exceeds this temperature point, the screen will display an error message "Er.H".	When you operate this device, the environment temperature cannot be greater than 40°C.
Er.L	This device has a lower limit of operating temperature of 10°C. When the ambient temperature exceeds this temperature point, the screen will display an error message "Er.L".	When you operate this device, the environment temperature can not be lower than 10°C.
Err	When the environment temperature changes rapidly about 5 degrees, it will show "Err" during measurement in object temperature mode, then shutdown automatically.	When it shows "Err", please place this product and keep steady in current environment for 30 minutes before next measurement.

 When the operation is wrong, **if the sound is turned on**, then the sound will read the value and remind at this moment. (A short beep will be heard).

TROUBLESHOOTING

Phenomenons	Reasons	Solutions
Device is powered on, but the screen does not light up.	The battery is exhausted.	Replace the new battery.
	The battery polarity is wrong.	The battery polarity is the same as the battery case.
The measurement temperature is low.	The measurement position is not correct.	Measure the temperature correctly according to the instructions.
	There are dirt stuff in the sensor or ear canal.	Please clear the dirt before the measurement.
Big temperature fluctuations with continuous measurement.	The measurement interval is too small.	The interval for each measurement should be above 10 seconds.

Product name:

Infrared thermometer for ear, forehead, object surfaces, and ambient.

Model number: Prontex X-Therm (REF: PG-IRT1603)

Product appearance dimensions: 31×175×72mm

Product weight: about 77g (except battery)

Measuring range:

- Ear: 34.0°C-43.0°C.(93.2°F-109.4°F)
- Forehead: 34.0°C-43.0°C.(93.2°F-109.4°F)
- Object temperature: 0°C- 93.2°C (32°F-199.7°F)

Resolution ratio: 0.1°C/°F

Accuracy: (35.0°C~42.0°C)±0.2°C
(95.0°F~107.6°F) ±0.4°F
Other temperature ±0.3°C.

Operating conditions:

- Temperature: 10.0°C~40.0°C(50.0°F~104.0°F)
- Relative maximum humidity: 15%RH~93%RH
- Atmospheric pressure: 70kPa~106kPa

Transportation/storage conditions:

- Temperature: -25°C~55°C (-13°F~131°F)
- Relative maximum humidity: 0%RH~93%RH
- Atmospheric pressure: 50kPa~106kPa

Display screen: LCD display screen, 4 bit numbers and special icons.

Sound:

when you turn on the device and ready to measure, a short beep will be heard. The measurement is finished with a long beep. System error or fault: short beeps for three times.

Fever alert:

short beeps for ten times come with urgency.

Memory: in memory mode, it can record nine temperature numbers.

Automatically shut down: if no operation for 30 seconds, it will shutdown automatically.

Battery: two pieces of 1.5V AAA batteries
(alkaline batteries are recommended to use).

Period of use: five years

Packing parts list:

1. Main body
2. Product manual ITA + ENG

Applied part: Type BF

Measuring Site: Forehead

Reference Site: Armpit

Measuring Site: Ear

Reference Site: Armpit

ADJUSTED MODE DECIMAL FIGURES

At the same time, press the “**ⓘ**” button + the “**⏸**” button for 5 seconds to display CH1, indicating that the ADJUSTED MODE has been entered. Press the MODE button again to display CH2, indicating that the accuracy is two decimal places.

Appendix 1 Guidance and Manufacturer Declaration Tables

Guidance and manufacturer's declaration—electromagnetic emissions		
Prontex X-Therm Infrared Thermometer is intended for use in the electromagnetic environment specified below. The customer or the user of Prontex X-Therm Infrared Thermometer should assure that it is used in such an environment.		
Emissions	Compliance	Electromagnetic environment-guidance
RF emissions CISPR 11	Group 1	Prontex X-Therm Infrared Thermometer uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	Prontex X-Therm Infrared Thermometer is used in home and it's powered by DC 3V.
Harmonic emissions IEC 61000-3-2	N.A.	
Voltage fluctuations/ flickeremissions IEC 61000-3-3	N.A.	

Guidance & Declaration — electromagnetic immunity			
Prontex X-Therm is intended for use in the electromagnetic environment specified below. The customer or the user of Prontex X-Therm should assure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment-guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
Radiated RF IEC 61000-4-3	10 V/m 80 MHz to 2.7 Ghz 385MHz-5785MHz Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communication equipment (Refer to table 9 of IEC 60601-1-2:2014+ A1:2020)	10 V/m 80 MHz to 2.7 Ghz 385MHz-5785MHz Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communication equipment (Refer to table 9 of IEC 60601-1-2:2014+ A1:2020)	Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of Prontex X-Therm, including cables specified by the manufacturer. Otherwise, degradation of the performance of this device could result.

WARRANTY

This device is guaranteed against manufacturing defects for a period of 2 years from the date of purchasing as shown on your receipt (maintain original copy to have right to warranty services). During this period, the unit will be repaired or replaced free of charge if the fault is due to faulty design or manufacture. This warranty does not cover damages or malfunctions caused by use improper or contrary to the instructions of this manual and is not applicable to the commercial batteries supplied at the time of purchasing.

CAUTION

Manufacturer/Distributor cannot be held liable for accidental or indirect damage, if changes have been made to the device, unauthorised repairs and/or technical interventions have been performed or any of its parts have been damaged by accident, improper use and/or abuse.

Any unauthorised intervention, even minimally, on the device voids the warranty immediately and, in any, case does not guarantee correspondence with the technical and safety requirements of the Dir. 93/42/CEE (MDD) and subsequent changes and related standards.

Numero Verde

848-580600

SPLIT CHARGE

Shenzhen Pango Medical Electronics Co., Ltd

Main Site: No. 25 1st Industry Zone, Fenghuang Road,
Xikeng Village, Henggang Town, Longgang District,
Shenzhen, Guangdong China

Additional site1: 2-4 Floor, No.5 Shanzhuang Rd.,
Xikeng Village, Henggang Town, Longgang District,
Shenzhen City, Guangdong Province, China

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Lotus NL B.V. - Address:
Safety S.p.A.